
CATALOG

Solutions to Laboratory Set-up by Meeting Today's Standard



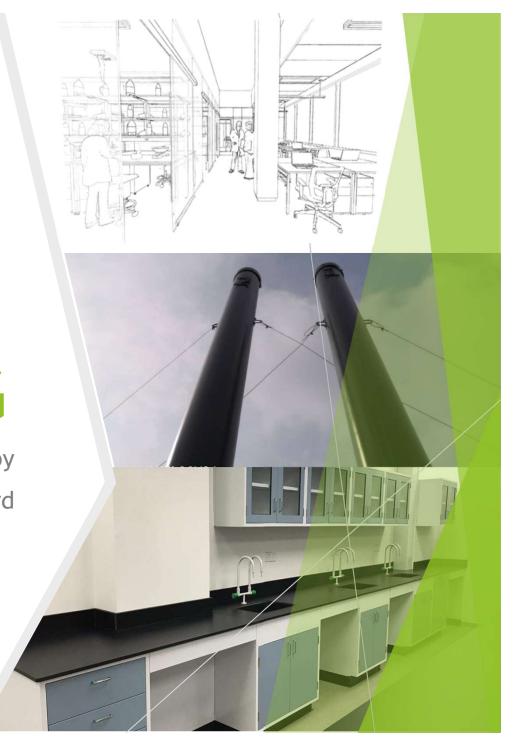

INTRODUCTION

Safemac Solutions Sdn Bhd is a Malaysian registered company that carries cutting edge technology scientific products. Established in 2010, being one of the business units under the **Safemac Group of Companies**, Safemac Solutions successfully developed its solid platform in various industries.

SAFEMAC SOLUTIONS is a leading one-stop laboratory design and facility setup provider.

SAFEMAC SOLUTIONS provides comprehensive laboratory facility design, renovation, laboratory furniture, products & speciality cabinets meeting today's demanding market of Laboratory solution.


Our design includes Chemical, Microbiology, Food & Beverages, Oil & Gas, Utilities, Biosafety Lab Grade, Clean Room facilities – Learn more from our design team.

Design & Installation Services for Local Exhaust Ventilation

We will design, supply, install and commission your LEV system from initial concept through to final completion with our experience within the industry. Our installation teams are experience in installing all types of industrial & laboratory ventilation equipment. LEV testing & Services are also available



Solutions for Laboratory Set-up

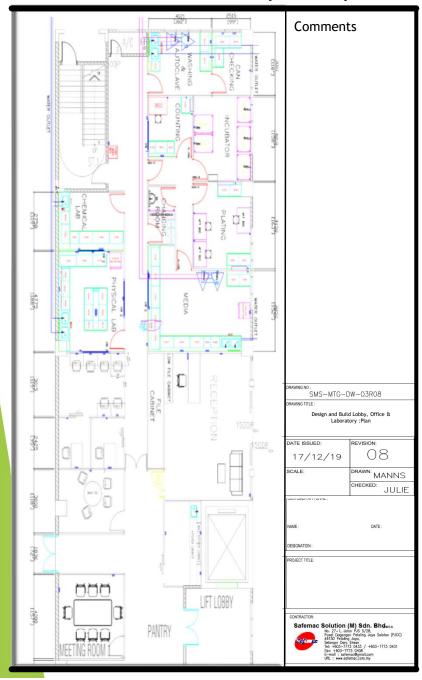













OUR PRODUCT RANGE

- **Laboratory Furniture** Page 5 - Pedestal - C Frame - Mobile - Custom - Clean Room Furniture
- **Local Exhaust Ventilation** Page 11 Fume cupboard **Extraction Arm System** Custom Made Exhaust Hood Testing & Services
- **Speciality Cabinets** Fume Cupboard Chemical Storage Cabinet Laminar Flow **Biosafety Cabinet**
- Cabin Laboratory Page 16 Gas Tank Cabinet

SAFEMAC SOLUTIONS Safemag

Solutions to Laboratory Set-up by Meeting Today's Standard

- **Safety Equipment** Page 21 Safety Shower
- Hand-held Eye-wash station
- Portable Eye wash Station
- Fire Blanket
- Fire Extinguisher
- Fire Alarm System
- Safety signage
- Page 24 Utilities
- Electrical bracket
- Service Riser
- Custom Sink
- Bin
- Laboratory Stool
- Reagent Rack
- Balance Table
- **Table Top** Page 29
- **Epoxy**
- Phenolic
- Stainless Steel
- Granite

Laboratory Furniture Standard & Custom



SAFEMAC SOLUTIONS











Steel Workstation Specification

Our metal cabinet are design and constructed from cold-rolled steel re-inforced, can be installed easily and can be re configured. All exposed cabinet surface are finished with epoxy coating for excellent chemical resistance.

Client can select a variety of selection from; Pedestal Free standing Modular system C Frame

Base cabinets are available in different colours of choice, doors, drawer fronts are factory installed to ensure product quality are maintain. Customer can also have the option to select type of handle for doors & drawer to suite their laboratory design.

Feastures:-

- Made of 1mm thick Electro-galvanise metal
- Removal Back Panel for easy services
- 1mm thick EG Epoxy EG Powder Coated Steel Body Structure
- Door flush with Super structure Frame
- Stainless 75mm 5 knuckle Hinges
- Roller Door Catch
- Door and Drawer Heads are "Double Pan" Construction and Sound Deadened
- Black Rubber skirting
- Heavi Duty Solid Steel Raceway 9rail) with nylon Tired Ball bearing Roller
- 2mm thick steel platform for adjustable leg mounting





Product Description

Wall Cabinet with 900mm height, 16 - 19mm thick Phenolic work surface with back splash, Standard width 750/900mm

Can have combination of sink, metal wall hung cabinets, knee space



Product Description

Island Cabinet with 900mm height, 16- 19mm thick phenolic work surface, standard 1200 - 1800mm width

Can have combination of sink, bin, reagent racks & knee space








Product Description - Modular

Wall Cabinet with 900mm height, 16 - 19mm thick Phenolic work surface with back splash, Standard width 750/900mm


Can have combination of sink, metal wall hung cabinets, knee space

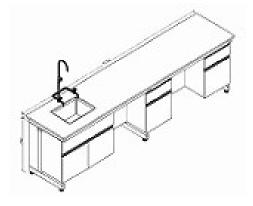

Product Description - Modular

Wall Cabinet with 900mm height, 16 - 19mm thick Phenolic work surface with back splash, Standard width 750/900mm

Can have combination of sink, Reagent Rack, knee space

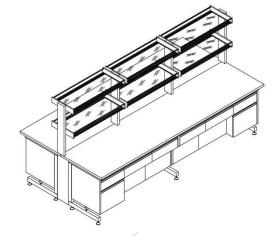



Mobile Cabinet



R&D Kitchen Lab



Custom Cabinet







Local Exhaust Ventilation Safemac Safemac SOLUTIONS

INTRODUCTION

Design & Installation Services for Local Exhaust Ventilation

We will design, supply, install and commission your LEV system from initial concept through to final completion with our experience within the industry. Our installation teams are experience in installing all types of industrial & laboratory ventilation equipment. LEV testing & Services are also available

Services Available

- Design & supply Extraction Hood
- Maintenance & Repair of Laboratory Hood
- Authority Application
- Analytical Instruments Hood supply & installation
- Supply Ductless hood for chemical cabinet

All our engineers & technician are train in house to assure customer satisfaction

























Speciality Cabinet Safemac SAFEMAC SOLUTIONS


OUR PRODUCT RANGE

- Speciality Cabinets
- General Purpose Fume Cupboard
- Ductless Fume cupboard
- Chemical Storage Cabinet
- Laminar Flow
- Biosafety Cabinet



Solutions to Laboratory Set-up by Meeting Today's Standard



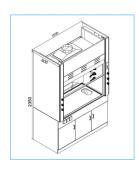

General Purpose Fume Cupboard

Model: 4ft (1219mm) Model: 5ft (1519mm) Model: 6ft (1829mm)



Laminar Flow

Size:
4 feet
5 feet

Model & Specification TBA



Ductless Fume Cupboard
4ft & 5Ft
Exhaust outlet with Prefilter
or Hepa Filter
2nd Stage activated carbon
Digital display
Airflow sensor
Fluorescent Light



Biosafety Cabinet

Size:
4 feet
5 feet

Model & Specification TBA




Walk In Fume Cupboard

Size: 6 feet

8 feet

Model & Specification TBA

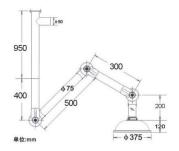

Chemical Fume Extractor



313 213 As 720 734

Custom Hood

Model & Specification TBA

All the fume cupboards are designed and build to cater different usage, application and environment. Focusing on high level of safety features and protection from hazardous substances, our fume cupboards are compliant with the standards such as EN14175 and Ashrae 110-95.

Our fume cupboards also come with uniform exhaust, air distribution in hood for safe and efficient handling of toxic, corrosive and other harmful fumes emitted from materials used in the laboratory.

Features:

-Double Wall Structure

Double wall structure consists of exterior housing (outer casing) constructed using baked epoxy steel panel which is easily detachable and chemical resistance to

external solvent and acid exposure. The internal hood (inner casing) uses one piece mold chemical resistance fiberglass for maximum space as well as easy cleaning.

-Air Baffle

Construction of fiberglass that is designed to maintain air flow of minimum point to point velocity variation with sash in maximum open position.

-Sash

Construction of tampering safety glass which provides clear and complete viewing of hood interior.

-Service Fixture

Fitted with service fixture for gases and water operated with the sash in any position and are controlled from outside the hood working chamber.

-Access Panel

Provided with access panels to facilitate easy removal for inspection or servicing of factory wired electrical components, control values and plumbing.

-Solvent Filter for Ductless Fume Cupboard

Hepa filter and carbon filter are used for dust and solvent fumes mounted on outlet duct.



Safety Equipment Safemac Safemac Solutions















Tall Case Glass Door - Metal

Size: 84-5/8" H x 22" D



Floor Standing
Fire Hose & Fire
Extinguisher Cabinet
- Metal



Tall Case Open Shelves - Metal

Size: 84-5/8" H x 22" D



Floor Standing
Gas Cylinder Cabinet
- Metal



Tall Case Metal Door - Metal

Size: 84-5/8" H x 22" D



Floor Standing Flammable Cabinet

- Metal Double Door
- Capacity 30, 40 & 60 Gallon





Floor Standing Acid Cabinet

- Metal
- Capacity 30, 40 & 60 Gallon

Accessories Safemae Safemae Solutions

Balance Table - Standard & Custom Size



Balance Table with granite top

- Standard Table height 900mm
- Anti-vibration purpose
- Top mounted shock & vibration absorbent insulation
- Comes with rubber base
 - Top Material Granite 40mm thickness (black/grey)
- Suitable for Analytical balance

Size: 36" x 24" X 36", 42" x 24" x 36"

Custom made available by enquiry only







Bin design







Service column design









Reagent Racks







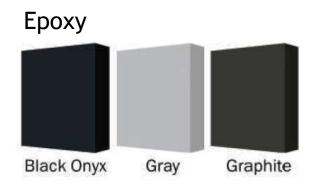






Table Tops Safemae SAFEMAC SOLUTIONS

Types of Table Top









Solicor-CR Features & Benefits



CHEMICAL RESISTANCE AESTHETIC APPEAL



Solicor-CR sets the industry standard for chemical resistance, providing designers and builders with a robust combination of performance and aesthetics.



Black core... No more. Gray compact with a gray core, and graphite compact with a graphite core set the new standard for an ultra-modern lab appearance.

FIRE RESISTANCE



Solicor-CR is rated CLASS A for fire and heat resistance, yet another attribute aligning it squarely with the laboratory

LIGHTWEIGHT



Solicor-CR worksurfaces are light enough to be recommended for both horizontal and vertical applications. Experience exceptional strength without all the weight.

MOISTURE RESISTANCE IMPACT & SCRATCH



Moisture resistance, even for the harshest acids, bases, solvents, and any other corrosive in a laboratory



Solicor-CR outperforms all competitors and sets a new standard for scratch and impact resistance. Full testing data available upon

COST-EFFECTIVE



Solicor-CR is cost advantageous, often requiring less upfront expense as well as lower life cycle costs. Guaranteed by our industry-best 10-year warranty.

USA MADE



Proudly manufactured & fabricated in the USA.

EASY TO INSTALL



At 7.4 lbs per square foot, the lightweight Solicor-CR is easy to transport, modify and install on the jobsite.

UV LIGHT RESISTANCE



exposure to direct sunlight, Solicor-CR coloration experienced unnoticeable variance from start to finish.




Using ultraviolet light to test prolonged

DURCON A WILSONART COMPANY


Worksurface Materials Comparison

Material	Chemical Resistance	Class A Fire Rating	Heat Resistance	Moisture Resistance	Impact & Scratch	UV Light Resistance	Edge Consistency
Epoxy Resin	~	~	~	~	~	~	~
Solicor-CR Lab Grad	e 🗸	~	~	~	~	~	~
Solid Phenolic Compa (SPC-CR) Lab Grade	ct 🗸	_	~	~	~	~	×
Competitor Lab Gra	de 🗸	~	~	~	~	×	_



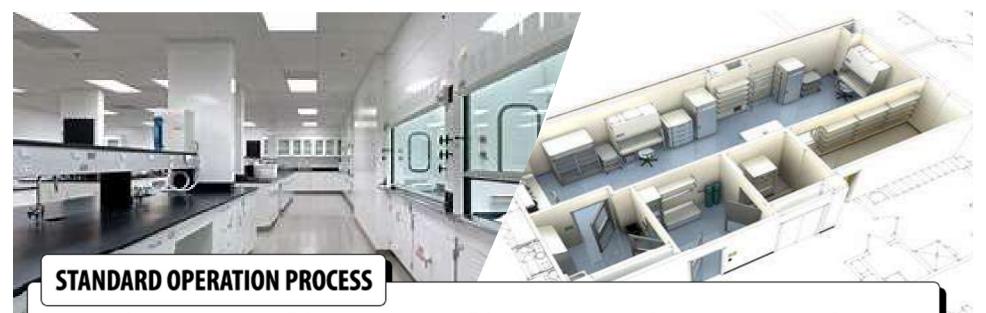



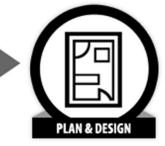











Safemac Solutions Sdn Bhd (898308-M)

No. 27-1, Jalan PJS 5/28

Pusat Dagangan Petaling Jaya Selatan 46150 Petaling Jaya, Selangor.

Tel: 03-773 0433 7773 0431

Fax: 03-7773 0438

Wedsite: safemacsolutions.com.my

Email:

sales@safemacsolutions.com.my